
Labofuge® 200

SERVICE MANUAL

P/N 12003630

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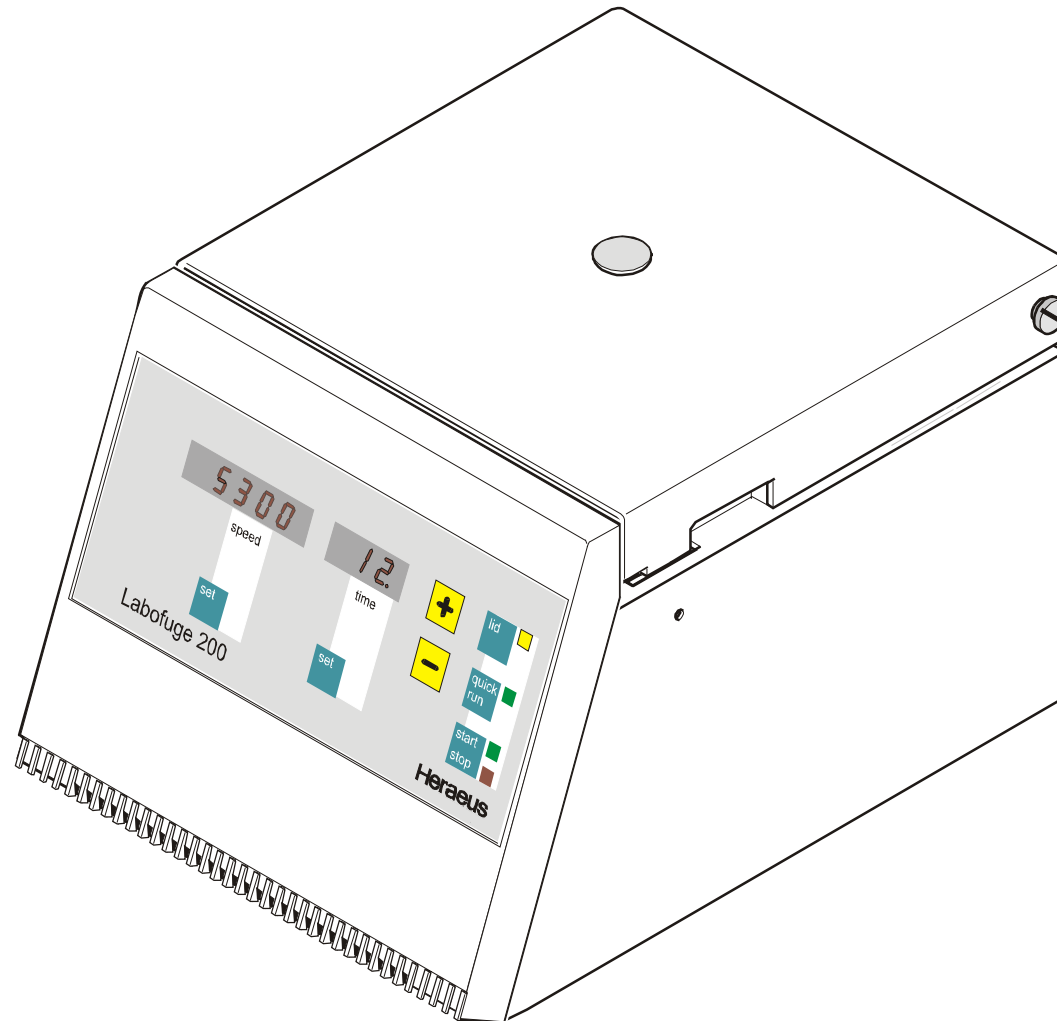TABLE OF CONTENTS

| Sect. | Title | Doc.- No. | Page |
|----------|---|-----------|---------|
| 1 | OPERATING INSTRUCTIONS (not part of this manual) | | |
| 1.1 | Brochure | Labo200_E | |
| 2 | SERVICE | Labo200_2 | |
| 2.1 | Servicing Schedule | " - " | 2-1/2 |
| 2.2 | Trouble Shooting | " - " | 2-3/5 |
| 2.3 | Test Points | " - " | 2-6 |
| 2.4 | Cleaning of Instrument Parts | " - " | 2-7 |
| 2.5 | Electrical Safety Check | " - " | 2-8 |
| 3 | FUNCTIONAL DESCRIPTION | Labo200_3 | |
| 3.1 | Block Functions | " - " | 3-1 |
| 3.2 | Functions of Main Board | " - " | 3-1/2 |
| 4 | DIAGRAMS | Labo200_4 | |
| 4.1 | Block Diagram | " - " | 4-1 |
| 4.2 | Wiring Diagrams | " - " | 4-2/3 |
| 4.3 | Wiring Connection Diagrams | " - " | 4-4 |
| 4.4 | Main Board 28 230V - Component Plan | " - " | 4-5 |
| 4.5 | Main Board 28 230V - Wiring Diagram | " - " | 4-6/8 |
| 4.6 | Main Board 178 120V - Component Plan | " - " | 4-9 |
| 4.7 | Main Board 178 120V - Wiring Diagram | " - " | 4-10/11 |
| 5 | DISASSEMBLY OF INSTRUMENT PARTS | Labo200_5 | |
| 5.1 | Housing / Casing Parts | " - " | 5-1 |
| 5.2 | Electrical Components | " - " | 5-2 |
| 5.3 | Drive Components | " - " | 5-3 |
| 5.4 | Mechanical Components | " - " | 5-4 |
| 5.5 | Service Kit | " - " | 5-5 |
| 6 | SPARE PART FIGURES AND LISTS | | |
| 6.1 | Break Down Drawings | Labo200_6 | 6-1/2 |
| 6.2 | Spare Part Lists | | |
| 7 | Preventive maintenance checklist | | |
| 8 | TECHNICAL BULLETINS | | |

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Printed in Germany – 3C 01/00 VN 4t Künzel ZF LF200 3d

LABOFUGE® 200 COMPACT CENTRIFUGE

Easy to Use

The Heraeus® Labofuge 200 is functional and sophisticated with an ergonomic design. It is ideal for use in medical practices, clinical and small laboratories and as a stand-by unit in large laboratories.

Benefits

- Easy to Use
- Able to run wide range of most popular medical tubes
- Maintenance free
- Safe

◀ Labofuge 200

Maintenance free

The Labofuge 200 is equipped with a brushless motor, so there is no need to replace carbon brushes. Beside saving time and money, this also ensures clean operating conditions.

Accessories

The Labofuge 200 comes complete with an autoclavable rotor. Made of impact resistant, fibre glass reinforced polyamide, this rotor offers high stability and outstanding run characteristics.

It is easily removed for cleaning in the laboratory washing machine. An extensive range of adapters permits centrifugation of all standard 5, 7, 10 and 15 ml tubes, including the popular "Monovette", "Vacutainer" and "VenoJect" blood collection tubes.

Safety


Samples processed in today's laboratories often harbor unknown risks. The Labofuge 200 complies with international safety standards and is equipped with a lid lock, lid interlock and steel armoured guard ring.

Functions

The microprocessor controlled Labofuge 200 is equipped with bright digital displays and touch-pad keys for quick and easy setting of speed and run-time. The values last used are stored. To repeat a run, simply press the start key.

It is also designed with a user friendly self diagnosis system which indicates faults directly on the display.


TECHNICAL DATA

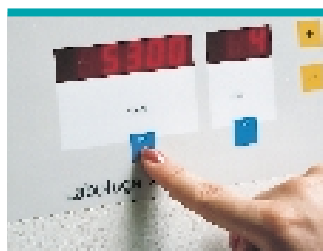Set the required speed...

...and set the required time

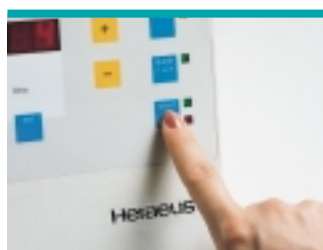With the "start" button, repeat runs can easily be recalled

| Labofuge 200 | 230 V version | 110-120 V version |
|---|--|---------------------|
| Description | Microprocessor controlled table top centrifuge | |
| Maximum speed (rpm) | 5300 | 5000 |
| Minimum speed (rpm) | 1600 | 1600 |
| Maximum RCF (x g) | 3030 | 2600 |
| Minimum RCF (x g) | 270 | 270 |
| Maximum capacity (ml) | 12 x 15 (glass) | 12 x 15 (glass) |
| Controller | Microprocessor controller for time and speed | |
| Drive | Brushless induction motor, microprocessor controlled | |
| Runtime | 1-99 min. and continuous operation (Hold) | |
| Program memory | Stores values last entered, unlimited | |
| Safety | Lid lock and lid interlock, steel armoured guard ring | |
| Design | Fibre glass reinforced polyamide housing with high noise and vibration insulation properties | |
| Dimensions (H x W x D) (mm / inch) | 240 x 284 x 375 / 9.5 x 11.2 x 14.8 | |
| Weight (incl. rotor) (kg / lbs) | Approx. 10.7 / 23.5 | Approx. 10.7 / 23.5 |
| Power consumption (W) | 65 | 65 |

ORDER NUMBERS

| Model | | Order No. |
|--------------|-----------------------|-----------|
| Labofuge 200 | 230 V; 50 / 60 Hz | 75003630 |
| | 110-120 V; 50 / 60 Hz | 75003631 |

| Accessories for Labofuge 200 | 230 V version | 110-120 V version |
|---------------------------------|---------------|-------------------|
| Angle rotor¹⁾ | 75003760 | 75003760 |
| Max. speed (rpm) | 5300 | 5000 |
| Max. RCF (x g) | 3030 | 2600 |
| Max. capacity (ml) | 12 x 15 | 12 x 15 |
| Max. radius (cm) | 9.65 | 9.65 |
| Acceleration time (sec) | 40 | 40 |
| Braking time (sec) | 45 | 45 |

Accessories for angle rotor 75003760

| Tube volume (ml) | Type of tube | Dimension (mm) | | Tubes per rotor | Type of adapter | Adapter Order No. |
|------------------|-------------------------------|----------------|--------|-----------------|----------------------|------------------------|
| | | Diam. | Length | | | |
| 7 | Glass tubes (DIN) | 12 | 100 | 12 | yellow | 75003227 |
| 15 | Glass tubes (DIN) | 16 | 100 | 12 | ⁴⁾ | - |
| 3-5 | VenoJect II (Terumo) | 13 | 75 | 12 | cream + rubber pad | 75003227 |
| 7 | VenoJect II (Terumo) | 13 | 100 | 12 | yellow | 75003227 |
| 9-10 | VenoJect II (Terumo) | 16 | 100 | 12 | sleeves | 75003763 ³⁾ |
| 3-5 | Vacutainer (BD) ²⁾ | 13 | 75 | 12 | cream + rubber pad | 75003227 |
| 7 | Vacutainer (BD) ²⁾ | 13 | 100 | 12 | yellow | 75003227 |
| 7 | Vacutainer (BD) | 16 | 75 | 12 | sleeves + rubber pad | 75003763 ³⁾ |
| 10 | Vacutainer (BD) ²⁾ | 16 | 100 | 12 | ⁴⁾ | - |
| 15 | Vacutainer (BD) | 16 | 125 | 6 | ⁴⁾ | - |
| 4 | Sarstedt Monovette | 11.5 | 83.5 | 12 | cream | 75003227 |
| 5 | Sarstedt Monovette | 13 | 90 | 12 | cream | 75003227 |
| 7.5 | Sarstedt Monovette | 15 | 92 | 12 | ⁴⁾ | - |
| 9 | Sarstedt Monovette | 16.5 | 92 | 6 | rubber pad | 75003762 |
| 12 (Urine) | Sarstedt Monovette | 16.5 | 101.5 | 6 | ⁴⁾ | - |
| Urine tubes | BD, Terumo | 16 | 100 | 12 | sleeves | 75003763 ³⁾ |

¹⁾ included with centrifuge

²⁾ tube incl. Hemoguard cap

³⁾ included with rotor

⁴⁾ not necessary

The adapter set 75003227 will contain the yellow and cream adapters and rubber pads.

For Ordering or Technical Information

Asia Pacific Kendro Laboratory Products (H.K.) Limited · Hong Kong · Tel. +852 2711-3910 · Fax +852 2711-3858 · info@kendro.com

Europe, Middle East, Africa

Kendro Laboratory Products International Sales · Hanau · Germany · Tel. +49 (6181) 35-300 · Fax +49 (6181) 35 59 44 · info@kendro.de

Latin America Kendro Laboratory Products International Sales · Newtown · USA · Tel. 1 (203) 270-2080 · Fax. 1 (203) 270-2210 · info@kendro.com

USA Kendro Laboratory Products · Newtown · Tel. 1 (800) 522-7746 · Fax. 1 (203) 270-2166 · info@kendro.com

Internet <http://www.kendro.com>

Quality Products – Lifetime Care

Registered to ISO 9001. Kendro Laboratory Products meet or exceed stringent quality and product safety standards: CE for the European Union, and UL, cUL or CSA standards for North America. ©2000 Kendro Laboratory Products. All Rights Reserved.

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2.1 Servicing Schedule (yearly procedure recommended)

2.1.1 Maintenance Routine without Dismantling the Centrifuge

2.1.1.1 Electrical Installation and Safety

- Switch OFF the centrifuge and disconnect the unit from power, check voltage supply and mains fusing (16 Amps, slow blow characteristic)
- Check condition of plug and wall socket - (let) replace defective parts
- Check cord condition and fixing / connection - replace or refit it

2.1.1.2 Location and Mechanical Installation

- Check the base (ground, table, lorry with lockable wheels etc.) For resonance-free and stable conditions
- Check for a well ventilated place and sufficient distances to walls or adjacent equipment, without exposure to direct sunlight
- Check the leveling of instrument (use a spirit level)

2.1.1.3 Lid Locking Mechanism and Safety

- Connect the centrifuge to power and switch ON
- Check for easy lid closing and self-acting lid opening - if in disorder, readjust lid's swivel hinge and/or lock assembly
- Check the central rubber gasket for lid sealing and replace it, if damaged
- For checking the safety: start the centrifuge, let it shortly run and stop it, the lid must not be unlocked by the microprocessor until the "end" message is shown on the display - if safety circuit is out of function, replace the main board

2.1.1.4 Cleanliness of Spin Chamber and Motor Casing

- Open the lid and remove the rotor (for loosening turn the central winged nut anti clockwise from the motor shaft)
- Clean the spin chamber with a dry and absorbent cloth (remove all dust and moisture - see also 2.5 Cleaning)
- Check the cleanliness of the motor casing and take care of the annular slot around the motor shaft: **penetrating fluids can damage the upper motor bearing**, remove fluids with an injector and/or absorbent paper

2.1.1.5 Rotor and Accessories Condition and Sealing

- Check the condition of rotors and accessory parts (especially all supporting or stressed partitions): the rotor and/or accessory parts must not be used any longer, if there are visible traces of mechanical damage

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Servicing Schedule**2.1.1.6 Rotor Fixing and Motor Shaft**

- Check the trouble-free condition of the locking nut and replace it in case of malfunction
- High performance rotors made of plastic have a limited service life. Due to safety reasons they have to be exchanged after 10,000 cycles or after being in use for five years. The expiry date can be seen on the date clock in the rotor lower shell. The two-digit number in the middle represents the year and the arrow points to the expiry month. The number of cycles may be checked at any time by closing the lid and keeping both “set” keys pressed. The accumulated rotor cycles will be displayed in the speed display.
- Check the condition of the drive motor shaft: the centrifuge must not be used any longer, if the drive shaft is damaged (bend, thread is worn out, horizontal grooves etc.)

2.1.1.7 Temperature Level

- Check the air inlet underneath the lid and under the bottom plate for free ventilation, insufficient air flow will lead to temperature rise of rotor, motor and electronic parts

2.1.2 Maintenance Routine after Dismantling the Centrifuge Casing**2.1.2.1 Motor Supporting Elements**

- Check the supporting and damping elements of the drive motor and replace them in case of increased rubber abrasion or abundance of imbalance but at least every 3 years

2.1.2.2 Braking Circuit

- Check the function of the braking circuit (even brake effect)

2.1.2.3 Lead and Screwing Connection

- Check the terminal and plug connections of all leads and on all boards and electrical components, tighten all loosen screwing connections, refit or replace defective parts
- Check the screwing connections of all boards, mechanical and electrical components and re-tighten them if necessary (use screw locking lac for motor mounts and lid lock assembly)

2.1.2.4 Protection Earth Core and Grounding Connections

- Check the protection earth core for continuity and all grounding plug connectors
- Check isolation resistance and accessible current (see 2.5)

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2.2 Trouble Shooting

| Error Indication | Error Cause | Possible Error Source | Corrective Procedure |
|---|---|--|--|
| Displays remain dark | No mains voltage supply | Mains fuse or circuit breaker failed | Check fuse or circuit breaker, replace or switch on again |
| | | Defective mains cord | Check instrument cord, replace defective parts |
| | | Defective unit fuse or fuses on main board | Replace it, if fuse blows again, disconnect electrical parts, otherwise replace main board |
| | No low voltage supply | Faulty connection from indication to power board | Check connections on CPU, indicat. Board and connecting leads, replace defective parts |
| | | Faulty indication or power board | Replace main board |
| Displays are illuminated, but drive doesn't start | Motor over-temperature switch has tripped | Motor temperature is higher than 120°C | Let motor cool down, then check temperature switch and leads with Ohmmeter |
| | Rotor didn't turn | Rotor is jammed | Check for easy rotor movement, remove any jamming objects |
| | | Motor is jammed | Replace motor |
| | Motor didn't start | Connections inter drive and main board | Check terminal and lead connections, replace faulty parts |
| | | Defective drive | Check resistance of motor windings, replace faulty parts |
| | | Faulty condenser | Replace condenser |
| | | Faulty main board | Remove main board completely and replace it |
| Drive makes noises-no good separation result | Mechanics | Wear out of motor rubber mount | Replace motor rubber mounts |
| | | Motor bearing | Change motor completely |
| | Electrical | Defective terminal connecton, faulty lead or motor winding | Check voltage on motor terminal and winding resistances -see test points on boards |
| | | Defective driving | Replace main board |

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Trouble Shooting

| Error Indication | Error Cause | Possible Error Source | Corrective Procedure |
|--|---|---|--|
| Drive doesn't decelerate | No brake current | Faulty main board | Remove main board completely and replace it |
| Lid cannot be opened by key pressure at standstill | Lid coil is not or not sufficiently supplied with voltage | Missing mains voltage | Remedy see above, manual opening only at standstill |
| | | PTC resistor has released | After a waiting time of 1-2 minutes press key again |
| | | Faulty driving or triac circuit | Replace the complete main board |
| | Faulty lid coil | Faulty winding of coil | Replace complete lid lock |
| Lid is not correctly locked | Lid bolt is jamming | Lid bolt is jamming | Push lid into lock and press the key again |
| | | Lid is deformed or disadjusted | Readjust the lid centrally |
| „Lid“ message appears in speed display | Lid was opened manually during run | Forbidden intervention Emergency opening can only be used at standstill | Close lid immediately, turn power off/on, wait for termination of br phase until end message appears |
| | Protection circuit (15V) for lid control was interrupted during the run | Defective micro switch or leads or connectors to micro switch are interrupted | Check leads and connectors to micro switch, in case of a faulty micro switch, replace lid lock device completely |
| „OPEn“ message in speed display | 15V supply circuit is interrupted at standstill | Defective micro switch or leads or connectors to micro switch are interrupted | Check leads and micro switch |
| „br “ message appears in speed display | Rotor comes to standstill without braking force | Short interruption of mains supply | Wait for rotor standstill (appr. 75 seconds) and re-start |

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Trouble Shooting

| Error Indication | Error Cause | Possible Error Source | Corrective Procedure |
|---|---|---------------------------------------|---|
| „E-“, „2“ message appears in speed field | Maintenance counter expired | Rotor has reached 10000 cycles | Replace rotor, NV-Ram and motor supports (Service Kit) |
| „E-“, „3“ message appears in speed field | No brake current | Faulty main board | Remove main board completely and replace it |
| „E-“, „12“ message appears in speed field | Checksum error of NV- RAM | NV-RAM is not initialized or false | Check NV-RAM and socket, insert the correct NV-RAM |
| | Disturbed data transfer from NV- RAM | Faulty main board | Remove main board completely and replace it |
| All dots in the display are illuminated | Maintenance counter expired | Rotor has reached 10000 cycles | Replace rotor, NV-Ram and motor supports (Service Kit) |

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2.3 Test Points

| Test Points | Unit value | Conditions |
|--|---|--|
| Mains terminal XA | 230V AC 120V AV | All given voltage and current values refer either to 230V ($\pm 10\%$) or 120V ($\pm 10\%$) |
| Lid micro switch Plug XB | 320V DC 160V DC | Voltage drop by open lid at 230V units Voltage drop by open lid at 120V units |
| Terminal XC Motor voltage | approx. 58V AC 105V AC 120V AC 80V AC 75V AC | Rotor #3760 not loaded, in each case measured after reaching the selected speed 1600rpm 4000rpm 5300rpm at the beginning of the braking phase at the end of the braking phase |
| Motor current I_M | approx. 1,2A 1,2A 0,8A 3,3A 2,0A | Rotor #3760 not loaded, in each case measured after reaching the selected speed with a soft iron or digital effective measuring instrument 1600rpm 4000rpm 5300rpm maximum during acceleration maximum during braking phase |
| Motor windings resistance 20°C -insulation value | 14,6Ω 14,6Ω > 10MΩ | switch OFF unit, pull off motor plug mains winding (black – green) auxiliary winding (black – yellow) resistance inter each phase and motor casing |
| Lid solenoid terminal XD | 105Ω 29Ω | switch OFF unit, resistance at 20°C (68°F) at 230V units resistance at 20°C (68°F) at 120V units |

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2.4 Cleaning of Instrument Parts

ATTENTION - WARNING!

The electrical and electronic components must not be cleaned with moist detergents!

For Cleaning the centrifuge housing or its accessories see Operating Instructions section 1 (maintenance and care).

- **Electronic components**

Clean dusty components carefully with a dry and soft brush and remove loose dust with a vacuum cleaner

- **Vent holes**

Remove dirt from vent holes of front panel or bottom plate using a brush and vacuum cleaner

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2.5 Electrical Safety Check

ATTENTION!

A final electrical safety check must be performed after each maintenance and/or repair!

- **Resistance check of protective conductor**


The measuring value of the resistance between the mains plug's grounding pin and the grounding conductors of the motor, electronic chassis and the casing must not exceed 200 mΩ.

- **Insulation resistance Check**

Check also the insulation resistance between the poles of the mains plug and the grounding conductor; the resistance value must be more than 2 MΩ.

- **Accessible current measured to EN 61 010**

The accessible current must not exceed 3.5 mAmps in single fault condition (interrupted protection earth wire)! In accordance with the EN61010, IEC1010 and UL3101 such a fault condition can be reproduced by the following measuring circuit.

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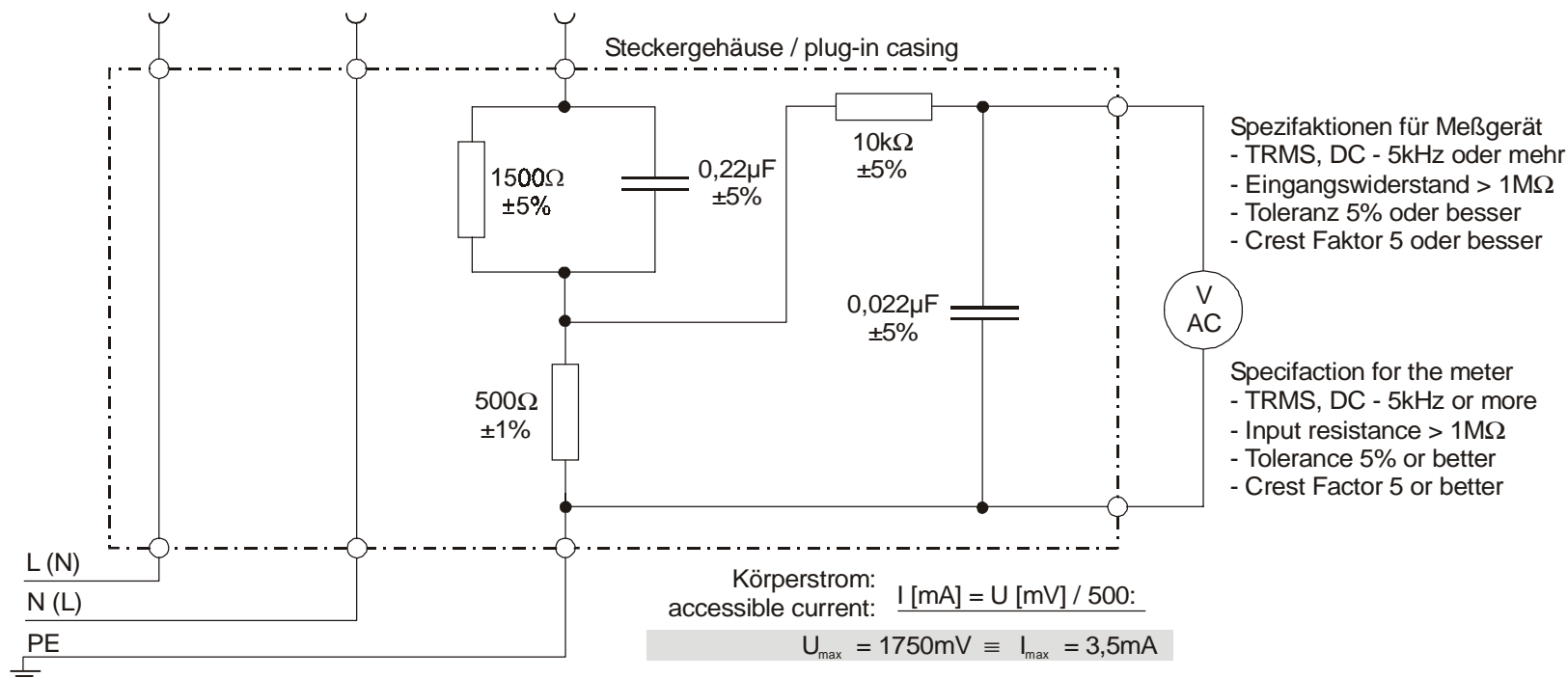3.1 Block Functions

The Labofuge 200 is a microprocessor controlled laboratory tabletop centrifuge with induction drive motor and integrated air cooling system.

The unit incorporates following boards and components (see block diagram 4-1):

- Main board 28 (for 230 V mains supply) or main board 178 (for 120 V mains supply) with microprocessor part and power electronics
- Key and indication board (programming: MEGACONTROL simple), part of the main board
- 1 phase induction motor with phase shift capacitor and integrated thermal over temperature switch (C. O. 125°C)
- Double lid lock assembly with solenoid and integrated micro switch (mechanical bolt keeper, magnetically dislocking), attached in front on both sides of the vessel

3.2 Main Board Functions

The main board 28 / 178 is mounted completely behind the front panel. The components on main board are arranged in following groups (see wiring diagram page 4-2).

- Board fusing (2 x 6.25 / 6,3 Amps. slow blow feature)
- Noise filter in accordance with EN 55011
- Power pack for low voltage supply of microprocessor part and power electronics (potential separated by transformer)
- Triac control circuit for lid solenoid (DC supply via bridge rectifier)
- Bridge rectifier for DC intermediate circuit supplying brake path and frequency converter
- Microprocessor part with CPU
- Exchangeable NV-RAM containing specific data of the unit

3.2.1 Power Pack

The power pack (transformer, bridge rectifier and voltage regulator) generates:

U1 = 5V: supplies central processor, key and indication board
reference potential: protective conductor (GND)

3.2.2 Intermediate Circuit with Brake Path and Frequency Converter

The DC intermediate circuit serves as an energy store between the AC power input and the transmitted motor performance. The intermediate circuit consists of a bridge rectifier (4 diodes) and two serial connected reservoir capacitors. The intermediate circuit is over voltage protected.

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Brake Path

To braking the centrifugal drive the induction motor is supplied with DC pulses from the intermediate circuit (affects as frequency zero). The pulse frequency is variable in dependence of the braking time (equal to speed), as a result of which the noise level is changed.

Frequency Converter


The frequency converter provides the drive motor with squarewaved voltage blocks having an amplitude of the half of the DC intermediate circuit voltage. These blocks are variable controlled in frequency and pulse-width modulation. Beginning with low frequency and small pulse-width blocks, both parameters will be increased up to the operators set speed. The set speed value will be reached with a tolerance of $\pm 5\%$ (within the valid mains voltage range)


3.2.3 Micro controller (Central Processing Unit)


The controller block includes the central processor unit (CPU) 87C51 and the data storage (NV-RAM) 9346. The ROM which contains the controlling program is integrated in the CPU. The most important operating parameters (e.g. the maximum speed or the last operator settings) are stored in the non volatile (NV) -RAM.


3.2.4 Key and Indication Board


The key and indication board is part of the main board. The connection to the main board is done by a 14 polar flat cable. The 7-segment displays, the control key LED's and the operating keys are managed by the CPU in multiplexing process. For the setting procedures serves MEGACONTROL "simple" (see Instruction Manual)

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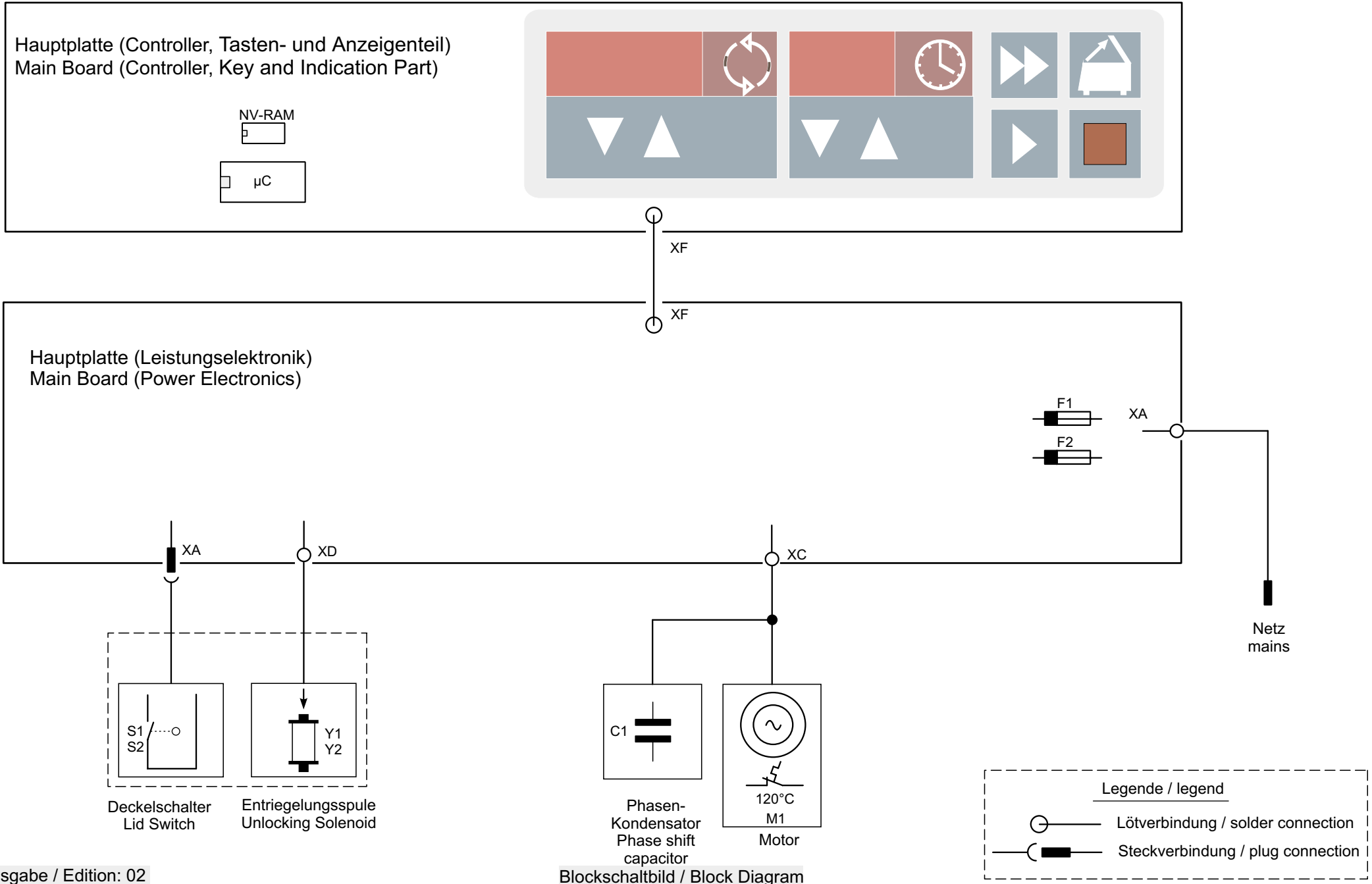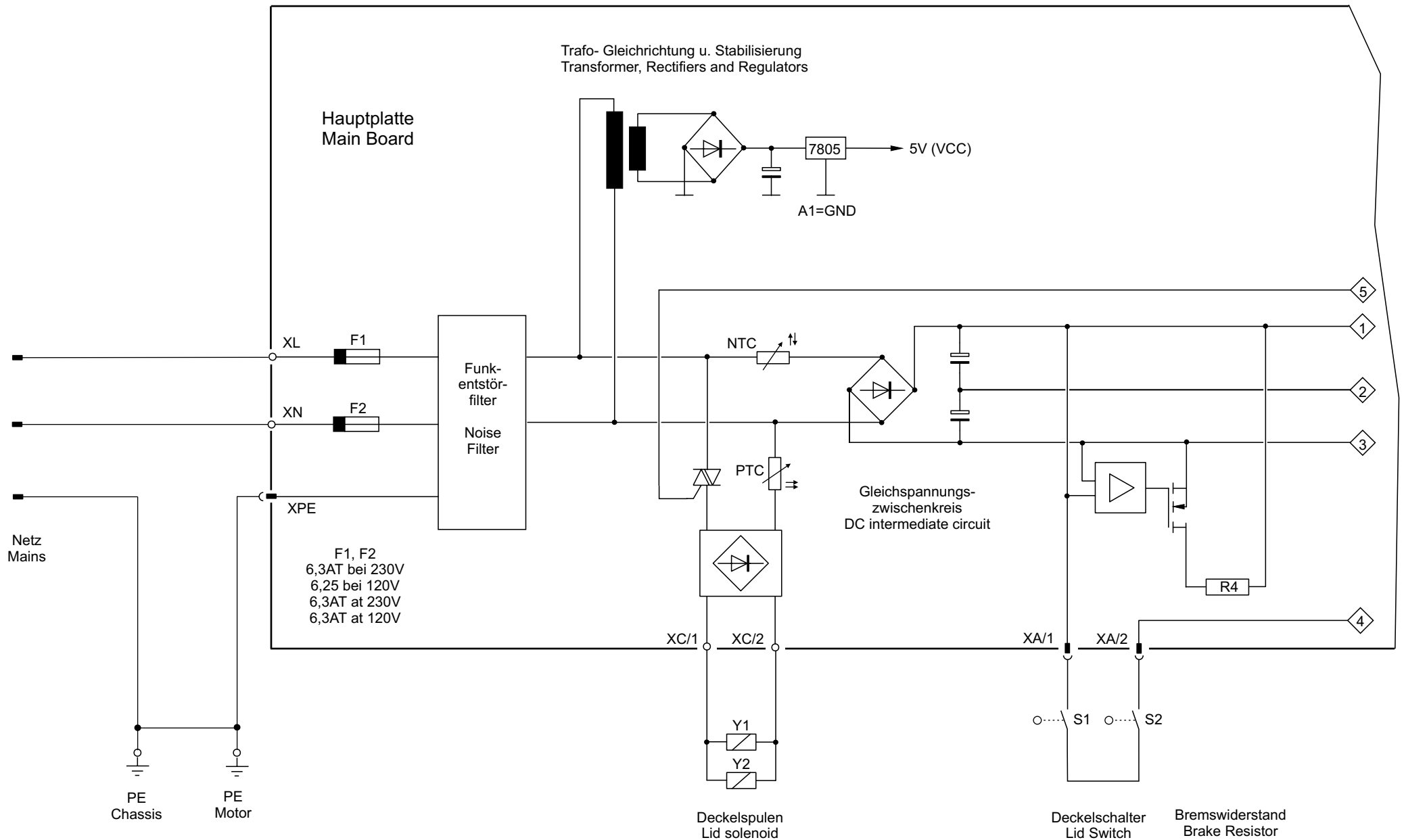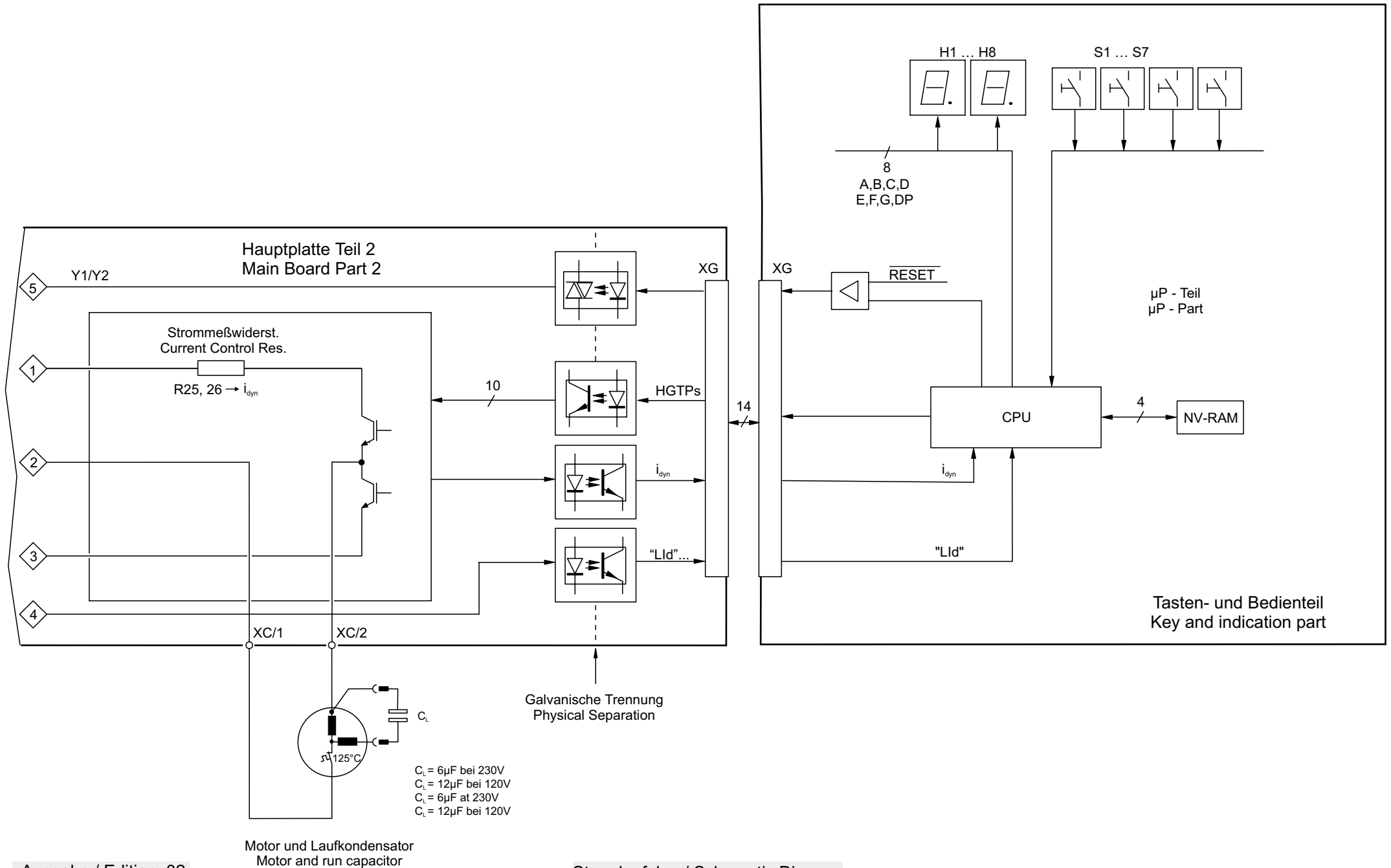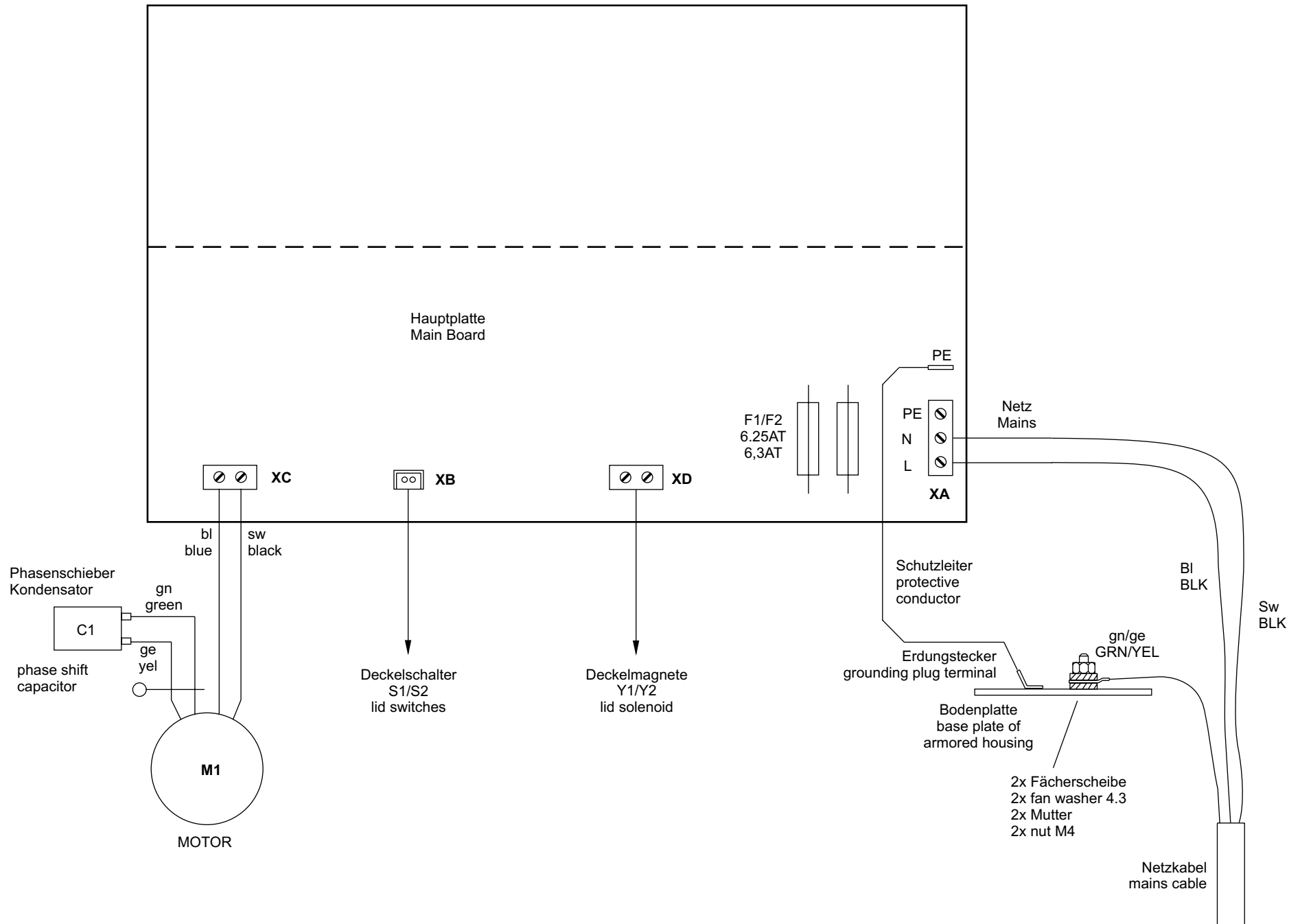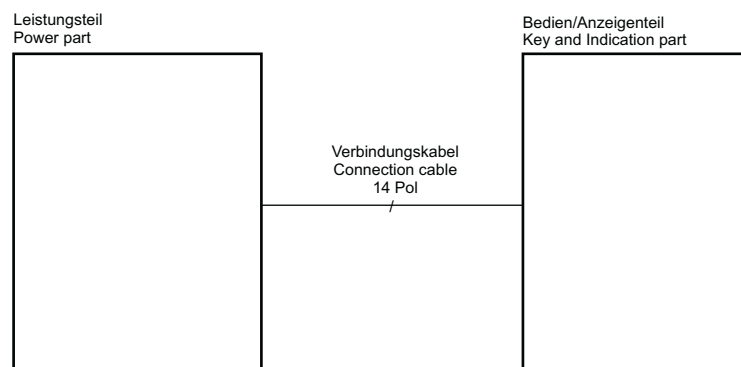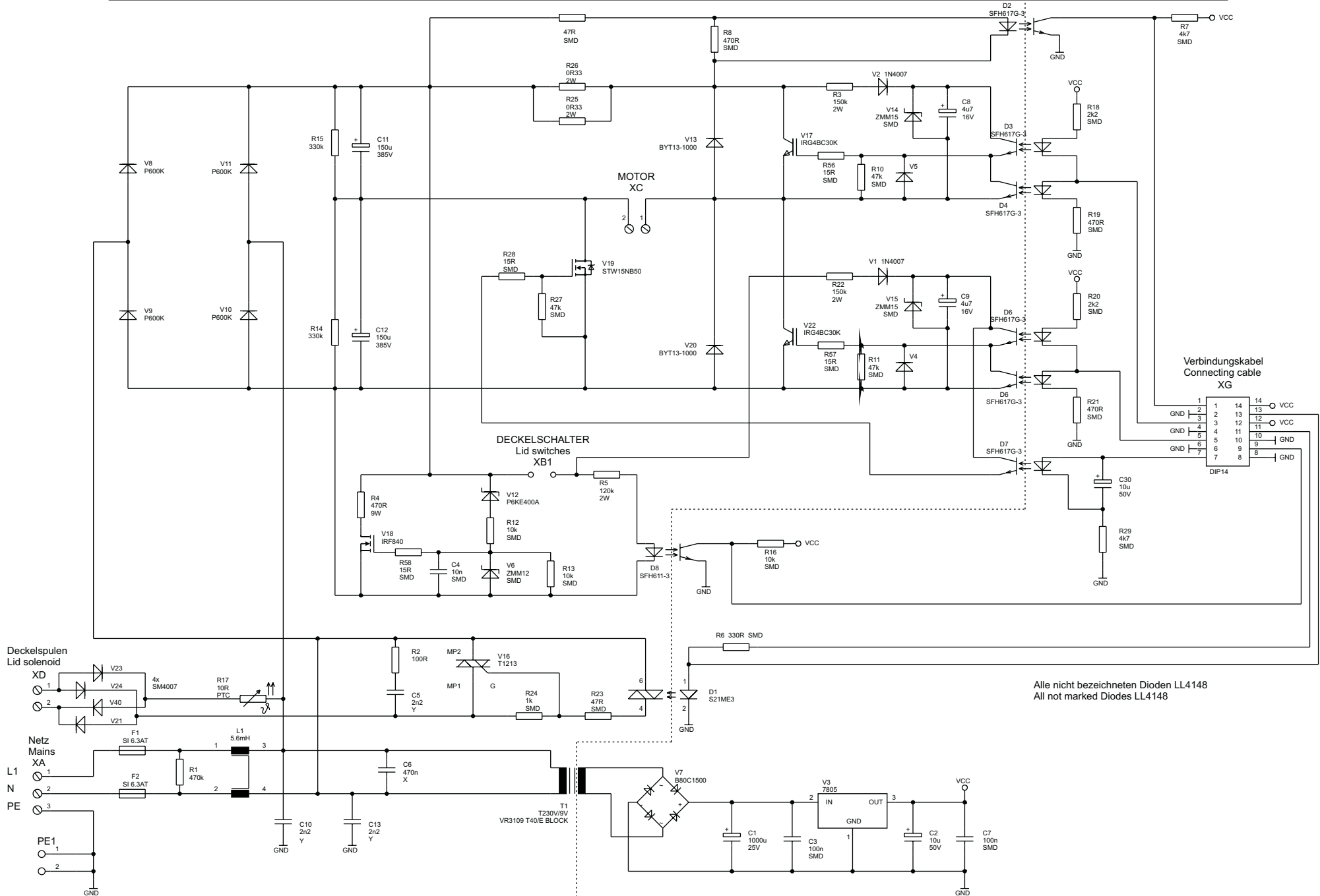
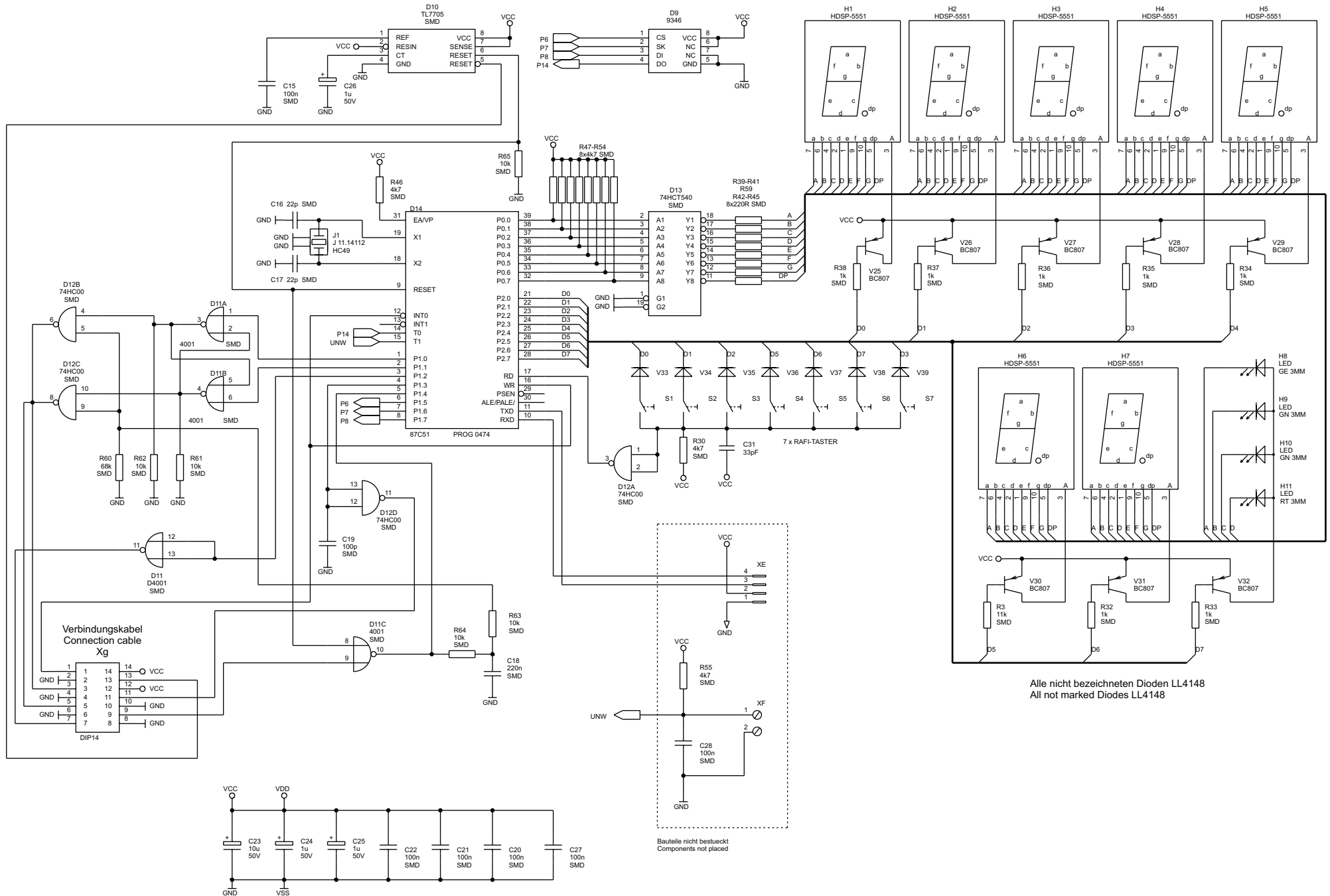
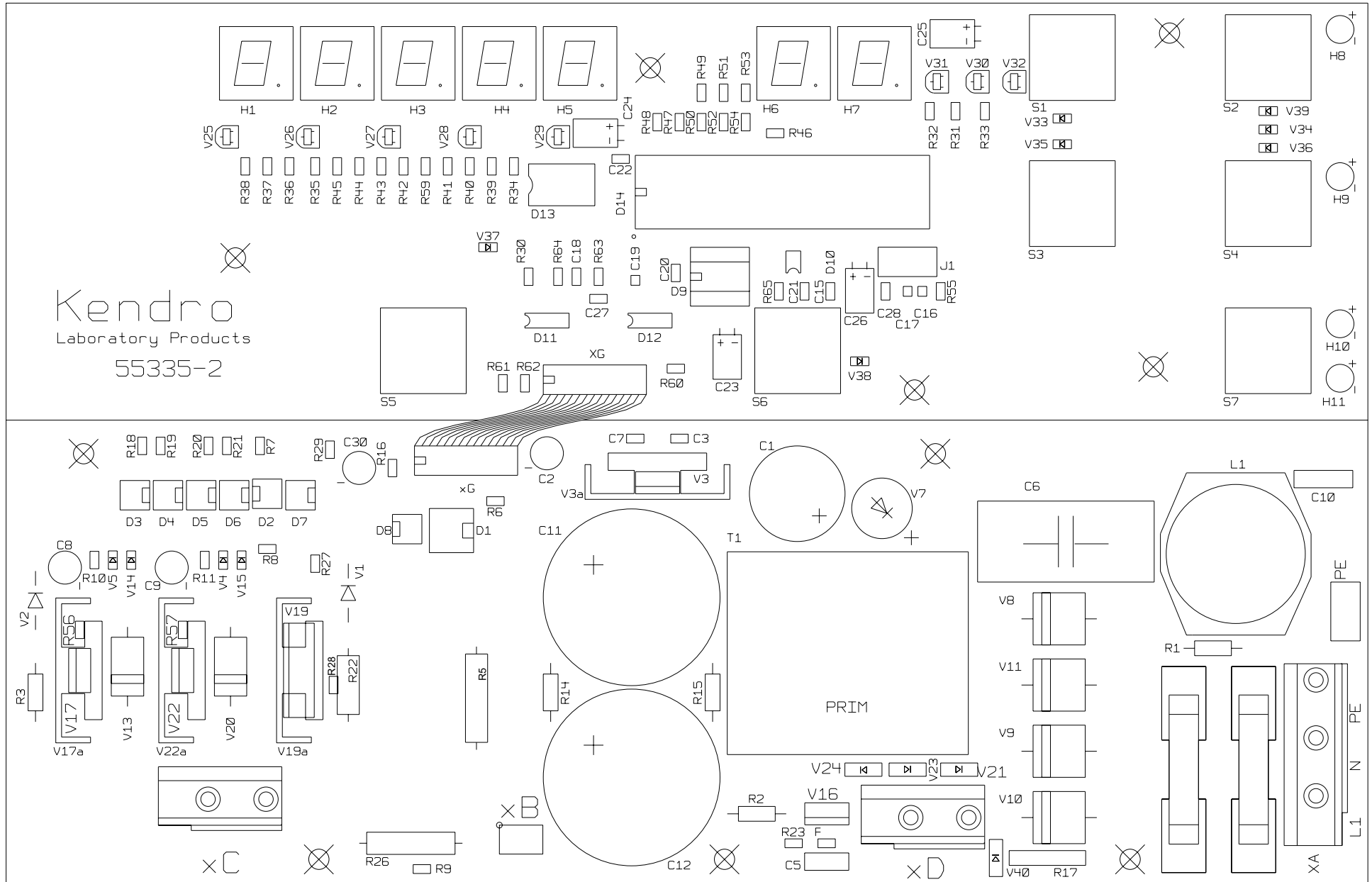
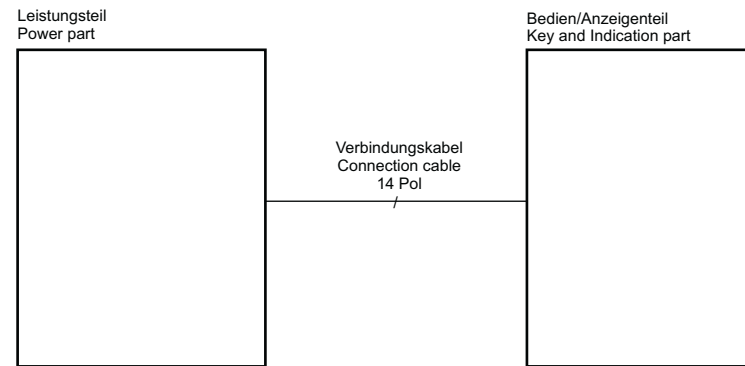
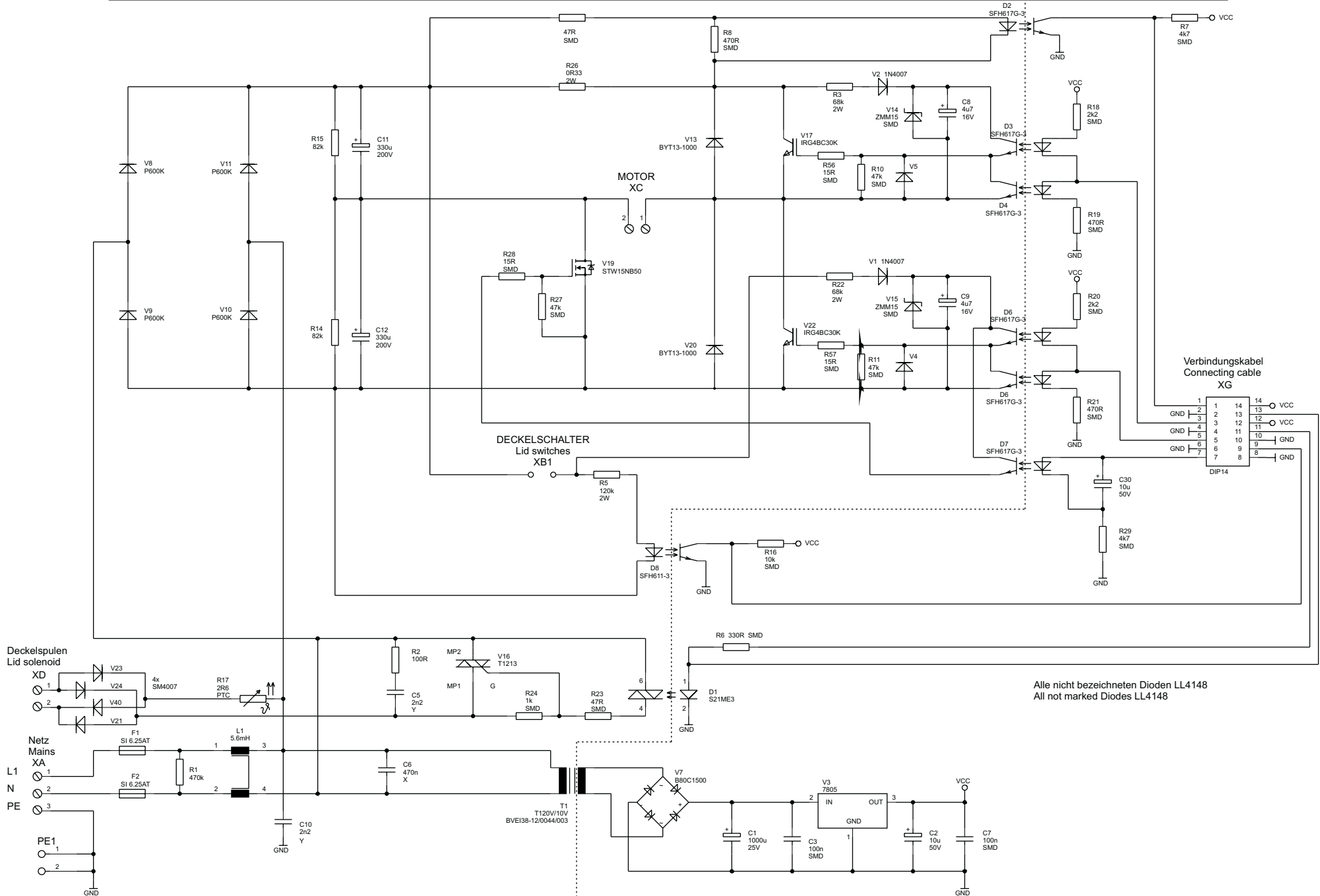
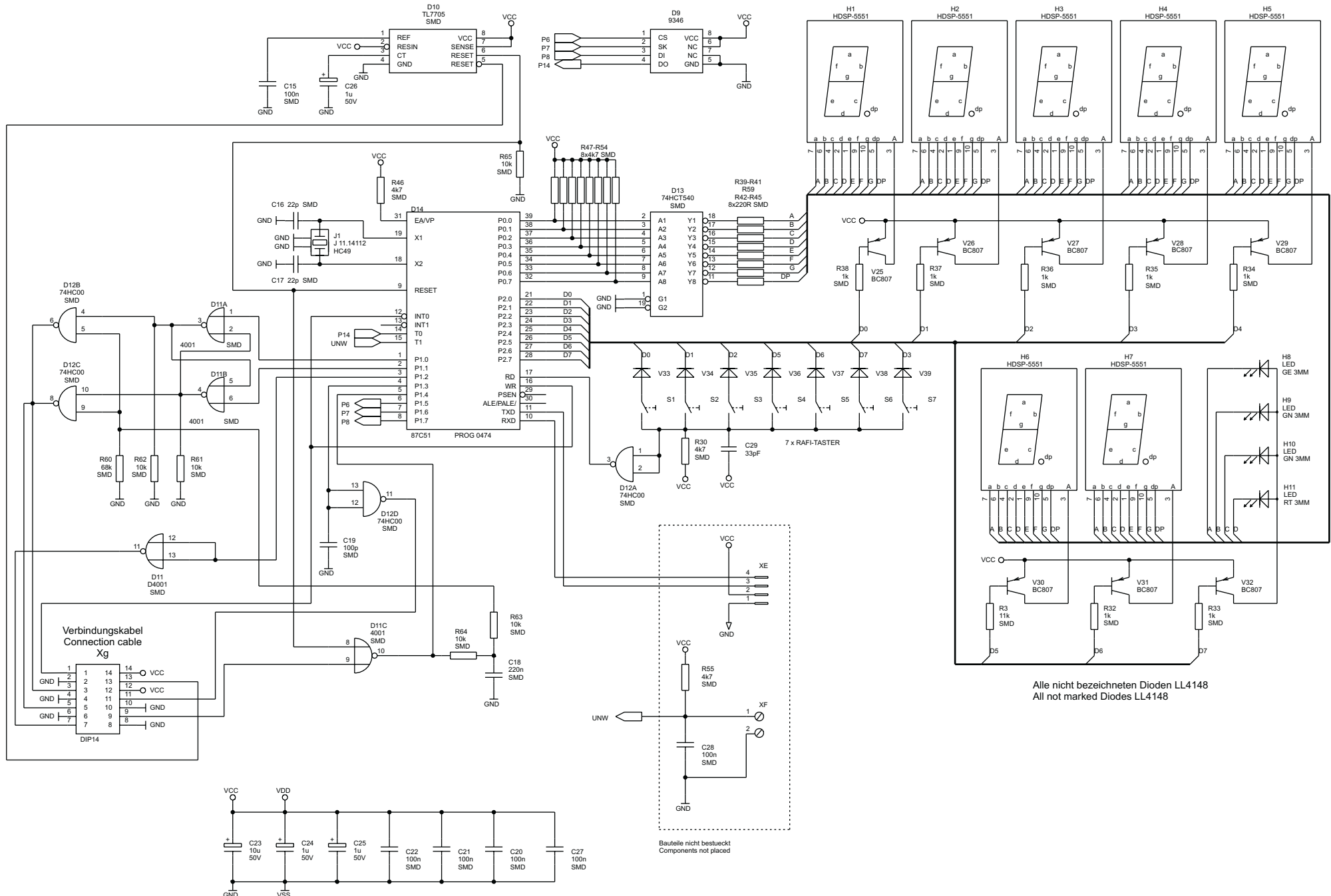
The index numbers stated in () reappear within the breakdown drawings and the spare part lists.

5.1 Dismantling the Housing

5.1.1 Front panel (102)

- Pull out the mains plug and place the instrument to the rear for the wide of the front panel
- Unscrew both screws (type Ejot 101) at the lower rim of front panel (use a small Phillips screw driver)
- Pull the panel's lower rim a little to front (but no more than 15° to avoid the braking off of the plastic hooks at the upper rim)
- Push the front panel in this position to the top to unlock it's casing connection and depose the panel in front of the unit
- Reassemble the front panel analogously in reverse order and press the lower rim tight to the casing when refitting the Ejot screws (but don't tighten the Ejot screws too much to avoid an overscrewing of the casing's plastic thread)

5.1.2 Lid (210)

- Open the lid and pull out the mains plug
- Unscrew both threaded hinge bolts (214), remove them together with bushings (212), spring washers (213) and rubber O-rings (215) and take off the lid
- Reassemble the lid analogously in reverse order (don't tighten the bolts too much, to avoid a squeezing off of the O-rings)

5.1.3 Dismantling the Casing (205)

- Open the lid and pull out the mains plug
- Unscrew the rotor and remove it from the motor shaft
- Dismantle front panel and lid - see 5.1.1 and 5.1.2
- Remove both Phillips screws (204) from bottom plate in front and lateral from armored chamber
- Disconnect the mains cable from terminal **XA** on main board and pull of the grounding plug
- Pick up the front panel, turn it to 90° and place it in front of armored chamber
- Press the casing upwards and remove it from the bottom plate passing the front panel in vertical position
- Reassemble the instrument analogously in reverse order

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5.2 Replacement of Electrical Components

5.2.1 Fuses on Main Board (113) replacement

- Dismantle the front panel (see 5.1.1)
- Check the micro fuses and if necessary replace them by equivalent and undamaged parts
- Perform a test run and search for blowing cause
- Reassemble the device analogously in reverse order

5.2.2 NV-RAM on Main Board (112)

- Dismantle the front panel (see 5.1.1)
- ATTENTION - CMOS components! Discharge your body before handling! Notice correct position of NV-RAM (8 pins) and pull carefully out of socket
- Reinsert the new NV-RAM correctly
- Reassemble the device analogously in reverse order reconnect it to power and perform a test run

5.2.3 Main Board Replacement (110)

ATTENTION! The key and indication board is part of the main board and cannot be ordered separately!

- Dismantle the front panel (see 5.1.1)
- Disconnect plug connections for lid switches **XB**, protecting earth conductor **PE**, and unscrew the lines for lid solenoid **XD**, motor **XC** and mains supply **XA**
- Remove all screwing of the main and indication board (Ejot type 111 and 122) with fiber washers (121) and take out both boards
- Touch a grounded receptacle to discharge your body before touching the sensitive CMOS components! Take a new NV-RAM (112) out of box (or re-use the old but trouble-free component) and insert it into the socket of the new main board (if necessary, remove a placed but non programmed NV-RAM before)
- Place the parts of the new board correctly and mount it with all screws and fibre washers
- Before assembling the front panel check the trigger function of all keys and correct the board mounting if necessary)
- Do not mix up disconnected cables during re-connection!
- Reassemble the device in reverse order and perform a test run, making sure the drive turns in the right direction (see imprinted arrow on rim of rotor chamber)!

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5.3 Replacement of Drive Components

5.3.1 Disassembly of Drive Motor (310)

- Remove the casing (see 5.1.3)
- Disconnect the motor leads from terminal **XC** and pull off the plug connectors from the capacitor (327)
- Take the leads out of the cable holder (308)
- Unscrew the 3 nuts (313) and remove them together with the lock washers (314)
- Lift the motor out of the rubber mounts and remove the motor from the bottom plate
- Reassemble the motor analogously in reverse order and tighten the 3 nuts with use of torque key (5Nm) and secure them with locking lac
- Reassemble the device in reverse order
- Perform a test run and check the sense of rotation (see imprinted arrow direction on casing)

5.3.2 Motor supports (316)

All 3 rubber supports (316) have to be replaced at the same time and at least every three years! By the way the 6 double-sided sandpaper discs (315) must be exchanged, too.

- Unscrew the drive motor (see 5.3.1), hold it with one hand and lay the armored chamber together with motor onto the side
- Pull the motor from the supports and lay it carefully onto the side
- Remove the 3 screws (224) from the bottom plate and take out the rubber supports together with sand paper discs (315)
- Install the new rubber mounts together with new sandpaper disks, put the lock washers in place, tighten the screws evenly and secure them with screw locking lac
- Put the unit again onto its feet, mount sandpaper disc, motor and lock washer and tighten the 3 nuts evenly by use of a torque key (5 Nm) and secure them with screw locking lac
- Reassemble the device in reverse order
- At last perform a test run

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5.4 Replacement of Mechanical Components

5.4.1 Lid Lock Assemblies (300)

ATTENTION!

Both lid lock assemblies inclusive micro switches are combined to one spare part and cannot be ordered separately!

- Dismantle the casing (see 5.1.3)
- Disconnect the leads of the lid solenoids from terminal **XD** and pull out the plug **XB** for the micro switches
- Unscrew both upper screws (305) of the lid lock assemblies' attachment from the inside of the armored chamber and the lower Ejoy screws (304) from the outer side
- Push down the levers of micro switches, guide the micro switches through the openings of the armored chamber and remove both lid lock assemblies completely with rubber base plates
- Reinstall the new lid lock assemblies analogously in reverse order and secure the thread of all upper fixing screws (305) with Loctite 221!
- Reassemble the device in reverse order
- Check the lid's correct closing and self-acting opening functions after the correct installation and replace the rubber profile (206)

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
5.5 Service Kit for Labofuge 200 / Medifuge 200


This service kit includes all parts which have to be replaced after termination of 10000 cycles. The kit includes different NV-RAMs which are destined for the use of the referring control program (signification of micro controller No.) on the specific main board.

ATTENTION! When changing the NV-RAM it is absolutely necessary to insert the right one into the main board with the right micro controller No.! Misuse can lead to main board destruction!

The specific variants can be identified as follows:

| Cat. No. of unit | Mains connection | NV-RAM signification | Controller signification | Additional identifying information and further notes |
|-------------------------------|--|-------------------------------------|--------------------------|--|
| 3630 3476 3627 15258 | 220V/50Hz 220V/50Hz 220V/50Hz 220V/60Hz | Labo 8 or Labo ALT | 0436 | minimum speed setting 3000 min-1, E-02 message, main board #150043 (components on both sides) |
| 3630 3476 | 220V/50Hz 220V/50Hz | 1052 V01 | 0465 | minimum speed setting 1600 min-1, no E-02 message, main board #150043 (components on both sides) or main board #150061 (SMD) |
| 3630 | 220V/50Hz | 1052 V03 | 0474 | cycle counter is set to 10500 the first time, main board 28 #150103 (SMD) or main board #150061 (older version V02 is cancelled) |
| 3631 | 120V/50Hz | 1066 V02 | 0474 | main board 14 #150076 or main board 178 #150178 |

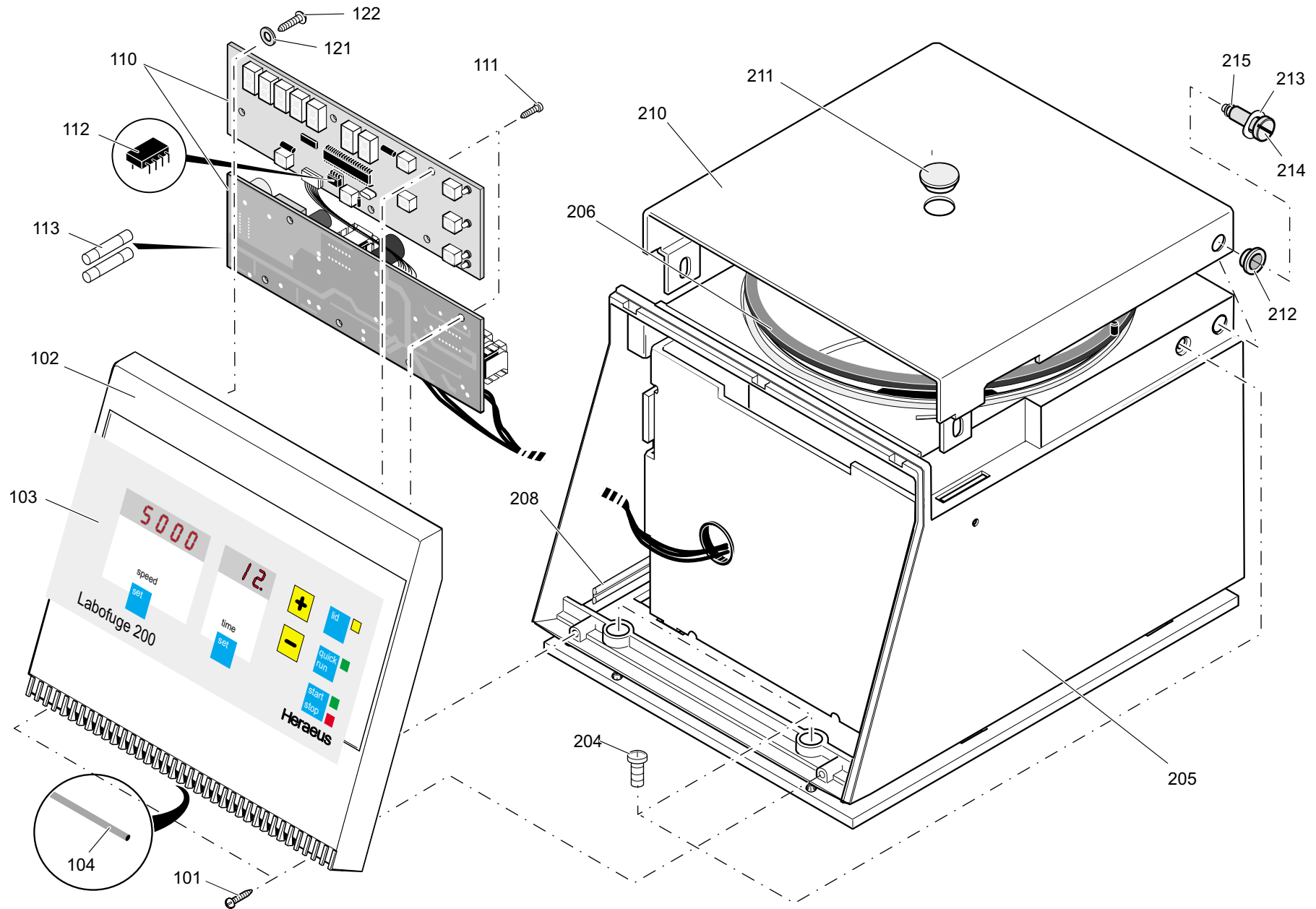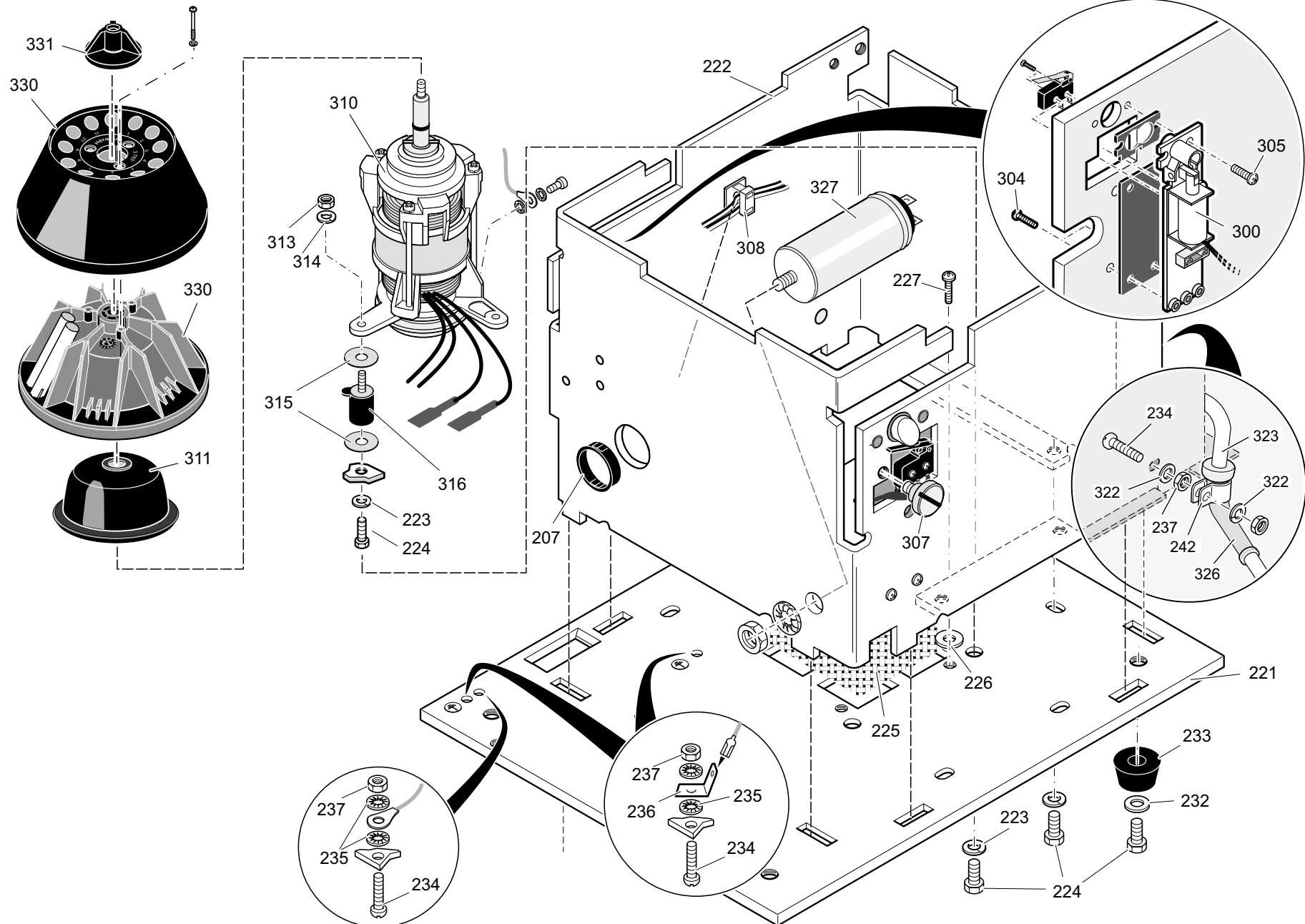75003630 _01

LABOFUGE 200 230V 50/60HZ

von Fabriknr.:

bis Fabriknr.:

| Index | Artikelnr. | Text |
|-------|------------|---------------------------------|
| 00101 | 20460133 | EJOT PT KB 030X014 SCHR LINS KR |
| 00102 | 20022598 | FRONTBLLENDE |
| 00103 | 20022604 | TASTERFOLIE LABO 200 HS |
| 00105 | 20290521 | MOOSGUMMI RUND 2,0 MM EPDM S |
| 00110 | 20150103 | HAUPTPLATTE 28 230V GS-VERSIO |
| 00111 | 20460137 | EJOT PT KB 030X008 SCHR LINS KR |
| 00112 | 70901052 | EE-PROM LABOF.200 220V 1052V03 |
| 00113 | 20230145 | SICHERUNG FEIN 6.3 A T (5X20MM) |
| 00121 | 20250031 | SCHEIBE 045X100X1,0 FIB |
| 00122 | 20460132 | EJOT PT KA 035X008 SCHR LINS KR |
| 00204 | 20460088 | D7985 050X008 SCHR LINS KR VA |
| 00205 | 20020726 | GEHAEUSE F. MEDIFUGE |
| 00206 | 20020717 | DICHTUNGSPROFIL F. MEDIFUGE |
| 00207 | 20220561 | SCHNAPPDURCHFHRUNG D19MM |
| 00208 | 20290473 | TESAMOLL 741 |
| 00210 | 70020819 | DECKEL RAL9002 BEDR.CAUT 2,5M |
| 00211 | 20050543 | SICHTFENSTER |
| 00212 | 20350211 | GLEITLAGERBUCHSE 8X10X5 FORU |
| 00213 | 20480250 | D0137 FEDERSCHEIBE A8 |
| 00214 | 20037565 | FLACHKOPFSCHRAUBE M.ANSATZ |
| 00215 | 20290419 | O-RING 0004,47X01,78 VITON (GU) |
| 00221 | 70022871 | BODENPLATTE 4 MM |
| 00222 | 70022870 | PANZERKESSEL GEBOGEN 4 MM |
| 00223 | 20480251 | FEDERRING DIN 128 A6 -1.4310 |
| 00224 | 20430216 | D0933 060X012 SCHR SKT VA |
| 00225 | 20022867 | DRAHTGEWEBE 0,5 MW0,8 60X60 (|
| 00226 | 20480224 | D9021 043 SCHEIBE VA |
| 00227 | 20460092 | D7985 040X006 SCHR LINS KR VA |
| 00232 | 20480157 | D0125 064 SCHEIBE VA |
| 00233 | 20049932 | FUSS |
| 00234 | 20460112 | D7985 040X016 SCHR LINS KR VA |
| 00235 | 20480108 | D6798 043 SCHEIBE F VA |
| 00236 | 20220510 | FLACHSTECKER 6.3X0.8MM 07F210 |
| 00237 | 20420076 | D0934 040 MUTTER SKT VA |
| 00300 | 20003620 | SCHLOESSER VOLLST. 230V |
| 00304 | 20460138 | EJOT PT KB 040X010 SCHR LINS KR |
| 00305 | 20440295 | D0965 040X016 SCHR SENK KR VA |
| 00307 | 20022873 | FLACHKOPFSCHRAUBE F.KLEINZE |
| 00308 | 20220416 | KABELCLIP |
| 00310 | 20210275 | MOTOR 211.30 220V/50HZ MEDIF. |
| 00311 | 20020738 | GUMMIFLANSCH F. MOTOR MEDI/BI |
| 00313 | 20420070 | D0934 060 MUTTER SKT VA |
| 00314 | 20480251 | FEDERRING DIN 128 A6 -1.4310 |
| 00315 | 10104080 | DOPPELS. SCHLEIFSCHEIBE 7X25 |

75003630 _01**LABOFUGE 200 230V 50/60HZ**

von Fabriknr.:

bis Fabriknr.:

| Index | Artikelnr. | Text |
|--------------|-------------------|--------------------------------|
| 00316 | 20022874 | GUMMIPUFFER RD 14X30, TAILLIER |
| 00322 | 20480130 | D0125 043 SCHEIBE VA |
| 00323 | 20190169 | NETZKABEL 3X1 MM EU |
| 00324 | 20220569 | BEFESTG.-SCHELLE F 9.5MM KABE |
| 00325 | 20480244 | FEDERRING DIN 128 A4 -1.4310 |
| 00326 | 20290189 | KNICKSCHUTZTUELLE HV 2107 (GU) |
| 00327 | 20540048 | KONDENS MAB M8 6/400 6UF |
| 00330 | 70901001 | SERVICE-KIT LABOFUGE 200 |
| 00331 | 20053358 | KNEBELMUTTER FUER 3760 |

Spare-Part-List

Kendro Laboratory Products

75003630_01

LABOFUGE 200, 230V +-10% 50/60HZ

Plant Osterode

from Serial-No.

to Serial-No.

| Index | Partno. | Text |
|-------|----------|-------------------------|
| 00101 | 20460133 | SCREW 3X14 |
| 00102 | 20022598 | FRONT PANEL |
| 00103 | 20022604 | KEY FOIL |
| 00105 | 20290521 | SEAL (FOAM RUBBER) |
| 00110 | 20150103 | MAINBOARD |
| 00111 | 20460137 | SCREW 3X8 |
| 00112 | 70901052 | NV-RAM LABOFUGE 200 NEW |
| 00113 | 20230145 | FUSE |
| 00121 | 20250031 | WASHER |
| 00122 | 20460132 | SREW 3,5 X 8 |
| 00204 | 20460088 | SCREW M5 X 8 |
| 00205 | 20020726 | HOUSING |
| 00206 | 20020717 | SEAL |
| 00207 | 20220561 | BUSHING |
| 00208 | 20290473 | SEALING |
| 00210 | 70020819 | LID |
| 00211 | 20050543 | WINDOW |
| 00212 | 20350211 | COLLAR |
| 00213 | 20480250 | SPRING WASHER |
| 00214 | 20037565 | SCREW |
| 00215 | 20290419 | O - RING 4,47 X 1,78 |
| 00221 | 70022871 | BOTTOM PLATE |
| 00222 | 70022870 | ARMED BOWL |
| 00223 | 20480251 | SECURING WASHER |
| 00224 | 20430216 | SCREW M 6 X 12 |
| 00225 | 20022867 | WIRE CLOTH |
| 00226 | 20480224 | WASHER A 4,3 |
| 00227 | 20460092 | ALLEN SCREW M4 X 6 |
| 00232 | 20480157 | WASHER |
| 00233 | 20049932 | FOOT |
| 00234 | 20460112 | SCREW M 4 X 16 |
| 00235 | 20480108 | FAN WASHER A 4,3 |
| 00236 | 20220510 | DOUBLE CONNECTOR MALE |
| 00237 | 20420076 | HEXAGONAL NUT M4 |
| 00300 | 20003620 | LID LOCK 230V |
| 00304 | 20460138 | SCREW 4X10 |
| 00305 | 20440295 | SCREW M4X16 |
| 00307 | 20022873 | SCREW |

Spare-Part-List**Kendro Laboratory Products**

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Kendro service Information

Heraeus Centrifuges

PREVENTIVE MAINTENANCE CHECKLIST

Labofuge 200
Medifuge 200

Type
S/N

Rotor Inspection

- Visual inspection
- Rotor locking screw
- Expiry date
.....
- Cycle counter checkup
.....

Preliminary Checks

- Cabinet
(clean, undamaged)
- Rotor chamber
(clean)
- Key foil
(clean, undamaged)
- Installation environment
(level & clearance)

Pre-run Checks

- Line voltage _____ V AC
- Display segments on power-up
- Function Keys
- Door interlock and latch
- Drive stop & Door unlatch (50 rpm)

Mechanical

- Motor bearings
- Antivibration mounts
- Motor rubber flange
- Chamber sealing

Performance Checks

- Acceleration time _____ seconds
- Instrument top speed _____ rpm
- Deceleration time _____ seconds
- Time (10 minute set) _____ actual min:sec

Safety test

- Accessible leakage current <3500 μ A
- Insulation resistance >10 Megaohm
- Earth Conductor Resistance <0.2 Ohm

Kendro service Information

Comments:

S.R. Signature _____
Date _____